GILSON® pipetman®

DIAMOND TIPS

PIPETMAN® DIAMOND TIPS

Standard is Gilson

Pipetman® Diamond Tips
Pure Performance Tips
for superl results

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HIGHEST QUALITY TIPS FOR SUPERIOR RESULTS

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MADE OF 100% PURE VIRGIN POLYPROPYLENE

No dyes. No additives. No trace metals. Smooth hydrophobic surface prevents liquid retention and ensures complete dispensing.

A PERFECT SEAL

Two solid sealing rings when used with PIPETMAN pipettes ensure more consistent and reliable results each and every time.

ERGONOMIC TIP FITTING

Best match between pipette and tips reduces the strain to fit and eject tips.

GRADUATED VOLUME MARKERS

Volume graduations make it easy to visually check the volume aspirated.

CONTAMINATION PREVENTION WITH PIPETMAN DIAMOND FILTER TIPS

Gilson's high-density polyethylene filter provides an effective barrier against aerosols.

FINE & FLEXIBLE POINT

The last droplet is easier to dispense.

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STERILIZED TIPS

- Sterilized with SAL 10⁻⁶ guaranteed in accordance with ISO 11137.
- Free of detectable RNases, DNases, DNA, RNA, ATP and pyrogens.
- All packs wrapped in medical device quality plastic that provides a complete shield against contamination.
- Sterilization expiration date on each rack.

Ordering Information

STANDARD PIPETMAN® DIAMOND TIPS

Sales unit \ Package								
		Racked		Racked (Sterile)		Bulk		Bulk
		TIPACK				ECOPACK		EASYPACK
	D10 0.1-10 µL • 30 mm	F171100	960 tips	F171101	960 tips	F161630	F161630—Bag 10*1,000 tips 1,000 tips	F161631 5*200 tips
	DL10 0.1-20 µL • 45 mm	F171200	960 tips	F171201	960 tips	F161450	F161450—Bag 10*1,000 tips 1,000 tips	F161451 5*200 tips
	D200 2-200 µL • 51 mm	F171300	960 tips	F171301	960 tips	F161930	F161930—Bag 10*1,000 tips 1,000 tips	F161931 5*200 tips
	D300 10-300 µL • 51 mm	F171400	960 tips	F171401	960 tips	F161730	F161730—Bag 10*1,000 tips 1,000 tips	F161731 5*200 tips
	D1000 100-1000 µL • 66 mm	F171500	960 tips	F171501	960 tips	F161670	F161670—Bag 10*1,000 tips 1,000 tips	F161671 5*200 tips
	D1200 50-1200 µL • 84 mm	F171600	960 tips	F171601	960 tips	F161110	F161110—Bag 10*1,000 tips 1,000 tips	F161111 5*200 tips
	D5000 0.5-5 mL • 10.4 cm	F161370 600 tips	F161370 -RACK 50 tips	-	-	-	-	F161571 3*334 tips
	D10mL 1-10 mL • 16 cm	F161230 240 tips	F161230 -RACK 30 tips	-	-	-	-	F161210 4*50 tips

EMPTY RELOAD RACKS For Eco & Easy Pack

DESCRIPTION	Reload Racks
D10	F171102 (2ea)
DL10	F171202 (2ea)
D200	F171302 (2ea)
D300	F171402 (2ea)
D1000	F171502 (2ea)
D5000	F161360 (1ea)
D10mL	F161220 (1ea)

EMPTY RELOAD BOXES For Refill Version & Reload Pack

DESCRIPTION
DL10/D200/D300 DFL10ST/DF30ST/DF100ST DF200ST/DF300ST
D1000/D1200/DF1000ST/ DF1200ST

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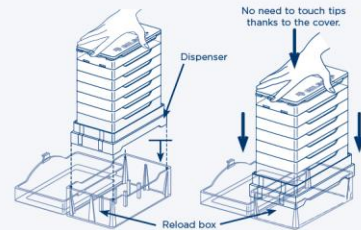
Pipetman® Diamond Filter Tips

DIAMOND FILTER TIPS Certificate of Conformity

The certification comprises the following tests:

Bio-molecule detection protocols

Disa-molecule detection protocols
For all tests, disposable tips are rinsed with sterile Dnase-RNase free water called «liquid extracts». Negative controls, positive and extract controls are made to validate all steps of the different assays.

Human DNA detection ≤ 2 pg

Human DNA detection < 2 pg
Preparation of sample with sterile and pyrogen free water. Extraction of DNA with a commercial kit. An aliquot of liquid extract is exposed to PCR reaction reagents containing primers (to amplify the β -actin housekeeping gene). A range of human DNA (2, 5, 10, 20 and 50 pg) is also amplified and used as positive control. Amplicons are obtained after cycles of amplification, analyzed by electrophoresis on a 2% agarose gel with TBE 0.5X buffer. No DNA has to be detected in the sample.

Received: 10 June 2009
Accepted: 17 July 2009

Bacterial DNA < 2 pg
Preparation of sample with sterile and pyrogen free water. Extraction of DNA with a commercial kit. An aliquot of liquid extract is exposed to PCR reaction reagents containing primer (to amplify rRNA 16S gene).
A range of bacterial DNA (2, 5, 10, 20 and 50 pg) is also amplified and used as positive control. Amplicons are obtained after cycles of amplification, analyzed by electrophoresis on a 2% agarose gel into TBE 0.5X buffer.

NO LHA HAVE TO BE d

RNA detection < 1 pg
Preparation of sample with sterile and pyrogen free water. Extraction, purification and RNA concentration are performed with a kit, and RNA converted to cDNA with an other kit, procedure recommended by the manufacturer. The cDNA is exposed to PCR reaction reagents containing primers (to amplify the β -actin housekeeping gene). A range of human RNA (1, 2, 4, 8 and 16 pg) is also amplified and used as positive control. Amplicons are obtained after cycles of amplification, analyzed by electrophoresis on a 2% agarose gel into TBE 0.5X buffer. No RNA have to be detected in the sample.

Endotoxin Detection ≤ 0.005 EU/ml

Endotoxin Detection < 0.005 EU/mL
Gibson sterilized DIAMOND Tips certified non-pyrogenic have been tested for bacterial endotoxins. Samples selected at random were tested and validated using the LAL kinetic chromogenic method D with < 0.005 EU/mL sensitivity. European Pharmacopoeia section 2.6.14 methodology for bacterial endotoxin testing.

DNoise detection ≤ 12.5 p.

DNAse detection (12.5 pg)
Preparation of sample with sterile and pyrogen free water. 10 μ L of liquid extract are incubated 2 hours at 37°C with variable concentrations of 100 plus DNA ladder (0, 10, 50 and 100 ng). Positive controls (0, 10, 50 and 100 ng of DNA + DNase I: 12.5 pg) and negative controls (0, 10, 50 and 100 ng of DNA) are incubated in the same conditions. Samples are analyzed by electrophoresis (2% agarose gel into TBE 0.5X buffer). The intensity of sample signals is compared with negative and positive controls. The degradation of DNA indicates the presence of DNase in the liquid extract. No degradation have to be observed, the degradation of negative control and the sample have to be identical.

RNAse detection < 0.25 µ

RNase detection < 0.25 ng
Preparation of sample with sterile and pyrogen free water: 10 μ l. of liquid extract are incubated 10' at room temperature (20-25°C) with 350 ng of RNA (0.1–20k RNA ladder). Positive control: (350 ng of RNA + 0.25 ng of RNase A) and negative control (350 ng of RNA) are incubated in the same conditions. Samples are analyzed by electrophoresis (2% agarose gel, 100 V, 100 mA, 1 h). The intensity of sample signals is compared with the positive and negative controls. The degradation of RNA indicates the presence of RNase in the liquid extract. No denaturation has to be observed. The intensity of control bands and the sample have to be identical.

ATB detection < 110% ms

ATP detection < 1.10⁻¹⁰ mol/L
Preparation of sample with sterile and ATP free water. Evaluation of ATP by luminescence (enzymatic reaction of luciferase - commercial kit). The intensity of the samples signals is compared with a negative control (DNase-RNase free water) and with ATP standards curve. The intensity of samples signals is compared with negative and positive standards curves. No ATP has to be detected.

Baseline distribution of 10 cases

Protein detection < 2 µg/mL. Preparation of sample with sterile and pyrogen free water. 150 µL of liquid extract are tested for the detection of proteins by absorption measurements (commercial kit). Positive control is a standard range supply by the kit, the negative control is sterile and pyrogen free water. Analysis are performed by absorbance reading at the plate reader. The intensity of samples signals is compared with negative and positive controls curve. 50 proteins have to be detected.

Inhibitors of PCR detection
Preparation of sample with sterile and pyrogen free water. Add of 50 pg of control genomic DNA for each sample. Amplification by PCR. Positive control is 50 ng of genomic DNA and an inhibitor standard (e.g. ethanol 0.25%, 0.50%, 1%, 2%, 5%), negative control is sterile and pyrogen free water. A DNA fragment with the same size of amplicon have to be detected.

This package contains a high-quality consumable manufactured in a controlled environment following ISO 14644-1. This lot is certified by validation of the manufacturing process following the requirements defined in ISO 18385.

The products are free of the following contaminants:

- Human DNA
- Bacterial DNA
- RNA
- ATP
- DNase
- RNase
- PCR inhibitors
- Endonucleases
- Proteases

This lot of tips has been Gamma irradiated for sterilization. The protocol analysis for validation of the sterilization process to ensure a Sterility Assurance Level (SAL) of 10^{-6} has been carried out following ISO 11173-2 for establishing the sterilization dosage.

Quality control and subsequent certification is done by independent, accredited laboratories (NF EN ISO/CEI 17025 - COFRAC).

Product certified ISO 18385
Company certified ISO 9001 and ISO 13485

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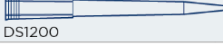
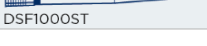
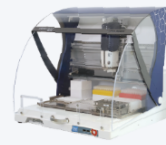
REFILL VERSIONS FOR PIPETMAX® CERTIFIED 384-WELL

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